

## **SPU High Back Slings**


#### **Features**

- Single Patient Use
- Full body support including head support and padded legs
- Use for bed, chair, commode & floor transfers
- Compatible with all market leading hoists
- Available with universal clip or loop attachments
- Safe Working Load:
  - S/M/L = 230kg
  - XL/XXL = 350kg
- Tested to Australian standards AS/NZS ISO 10535:2011
- Comprehensive label includes:
  - Unique serial number
  - QR code to download user manual
  - Space for patient ID
- Wetness indicator in commode region
- 'DO NOT USE' indicator exposed by laundering or immersion
- Clear size indication via label and universal colour-coded binding
- Unique rainbow webbing clearly differentiates from washable slings
- Heavy duty stitch-bonded fabric features include:
  - patient comfort soft polyester feels like clothing
  - breathability protects patient skin integrity
  - absorbency reduces fluid run-off
  - rip-stop nylon mesh grid for additional security
  - able to be 'spot cleaned' and spray disinfected
  - can be used in wet areas without discarding


## **NSW CF820 State Contract Approved Products**


### **SPU Clip Slings**



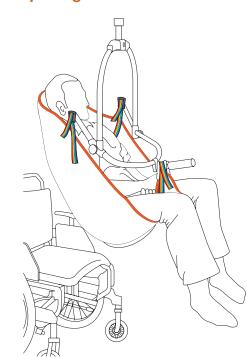

#### **SPU Loop Slings**

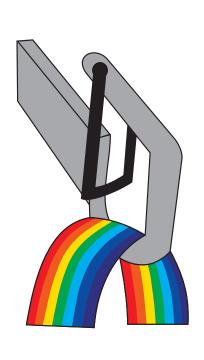


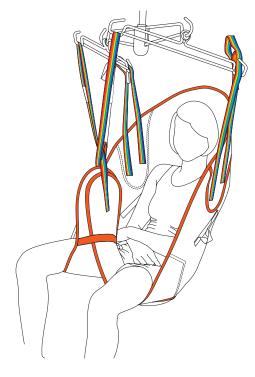



# **SPU High Back Slings**

# **GENERAL GUIDANCE - SLING FITTING & APPLICATION Sling Size Selection - Sling Dimensions**

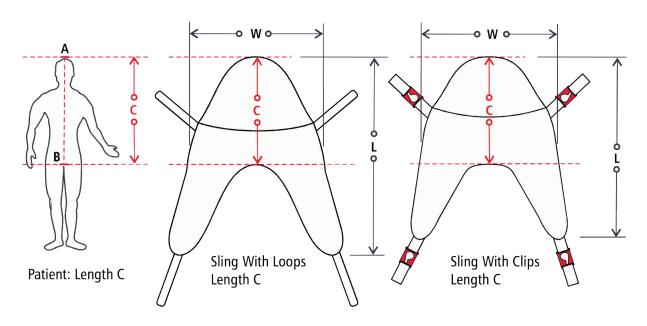

High Back Slings With Head Support: With Loops - Dimensions in CM Tolerence: ± 5%										
Product Code	Size	Sling Width W	Sling Length L	Sling Length C	Binding Colour	Safe Working Load	Box Quantity			
RX-EASLG-10-FB-LDS-S	S	78	109	75	Red	230kg	10			
RX-EASLG-10-FB-LDS-M	М	93	128	85	Yellow	230kg	10			
RX-EASLG-10-FB-LDS-L	L	111	153	97	Green	230kg	10			
RX-EASLG-10-FB-LDS-XL	XL	130	178	110	Blue	350kg	5			
RX-EASLG-10-FB-LDS-XXL	XXL	143	204	123	Orange	350kg	5			

High Back Slings With Head Support: With Clips - Dimensions in CM Tolerence: ± 5%										
Product Code	Size	Sling Width W	Sling Length L	Sling Length C	Binding Colour	Safe Working Load	Box Quantity			
RX-EASLG-10-FB-KC2-S	S	86	116	80	Red	230kg	10			
RX-EASLG-10-FB-KC2-M	М	99	133	89	Yellow	230kg	10			
RX-EASLG-10-FB-KC2-L	L	113	147	102	Green	230kg	10			
RX-EASLG-10-FB-KC2-XL	XL	123	173	118	Blue	350kg	5			
RX-EASLG-10-FB-KC2-XXL	XXL	137	192	130	Orange	350kg	5			

**Distributing for Avant Innovations** 

T: 1300 106 003 Dir: 07 3117 9539 E: sales@hmsedical.com.au

