




Experience true freedom and infinite mobility with the Magnum LED Surgical Headlight

IsoLux Magnum surgical headlamp gives total control to the surgeon.

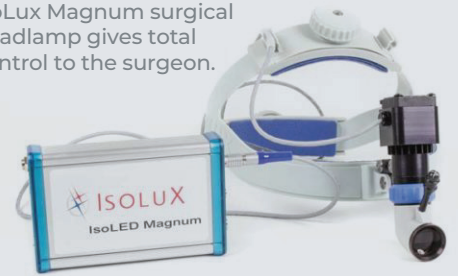

IsoLux Magnum Storage Case



# Isolux Magnum LED Headlight

## The Best LED Battery Powered Surgical Headlight for every situation

What makes the Isolux Magnum LED Battery Powered Surgical Headlight the best surgical headlamp system is the ability to adapt to the needs of Neuro, Plastic, Orthopedic, and Thoracic Surgeons. Magnum has been meticulously designed to offer reduced weight, optimal brightness, and color. The blue aperture control on the LED head gear allows for spot size adjustment during medical procedures.

The Magnum battery powered surgical headlight system is a feather weight solution that reduces neck strain caused by heavier headlamp systems. At only 10oz the Magnum is one of the lightest surgical headlight systems available. Trial a Magnum headlamp today to experience the difference.

Three LED intensity settings are available on the Magnum 30K Lux, 50K Lux, and 75K Lux. Each Lithium battery can last 9 hours on low, 6 hours on medium, and 3.2 hours on high intensity. Since the Magnum surgical headlamp comes with two Lithium batteries, the total operating time can range from 7 hours to 18 hours.

### Isolux Magnum Surgical Headlight Features

- 10mm - 100mm Variable Spot Control
- 10 oz Ultra Comfort Headlight Assembly
- Whisper Quiet LED Cooling Fan to Extend Runtime
- 5000K True Color Temperature
- 30,000 Lux, 50,000 Lux, & 75,000 Lux Intensity Settings

## HEADLIGHT SPECIFICATIONS

<b>CRI (Color Rendition Index)</b>	70
<b>Battery and Holster Weight</b>	12oz / 350 grams
<b>Color Temperature</b>	5000°K
<b>Headlamp Weight</b>	10oz / 300 grams
<b>Power Pack Operating Time</b>	3.2hrs @ High
<b>Three Intensity Levels</b>	30K Lux, 50K Lux & 70K Lux, all @12" / 30cm
<b>Power Packs Per System</b>	2
<b>Light Output @ 12"</b>	75K Lux
<b>Joy Stick</b>	Detachable for easy spot movement
<b>Quick Recharge Time</b>	2 hours, for 8 hours of maximum intensity
<b>Operating Time</b>	7 hours or continuous operation at Hi intensity
<b>Headband</b>	Lightweight Comfortable with washable foam inserts
<b>Warranty</b>	5 Years on LED Headlight 1 Year on Battery Pack